VEVOR

Upgrade · The Home Creator Way

Manual Rolling Walker

MODEL: HFK-9202E

VEVOR
Upgrade·The Home Creator Way

Manual Rolling Walker

MODEL: HFK-9202E

Note:The product picture is for reference, the actual details shall prevail

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

	Please retain instructions for future reference Protect yourself and others by observing all safety information, warnings and cautions. Failure to comply with instructions could result in personal injury and/or damage to product property.
	WEEE symbol. Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your Local Authority or local store for recycling advice.

PRODUCT STRUCTURE DIAGRAM

1.Armrest set
2.Backrest
3.Seat cushion
4.Hanging bag
5.Cup holder
6.Cane box
7.10-inch PU front wheel
8.8-inch PU rear wheel
9.Main frame set

APPLICATION SCOPE

1.Indoor use: This four-wheeled vehicle is an essential item for individuals with limited mobility to move freely at home.
2.Outdoor use: Effectively meets the needs of individuals with limited mobility to breathe fresh air outdoors or travel short distances.
3.Meets the user's need for short breaks.

MAIN PARAMETER

Model	HFK-9298
Wheel	8-inch PU Wheel
Load Capacity	300 lbs
Width	26inch

ASSEMBLY

Pre-Assembly Preparation

After unpacking the carton, first use scissors to remove the packaging. Please take care to protect all the parts inside the carton. If assembling on the floor, lay down kraft paper or a carpet first to prevent scratches. Then remove the inner carton and take out the parts. The carton contains loose components, as shown in Figure 1.

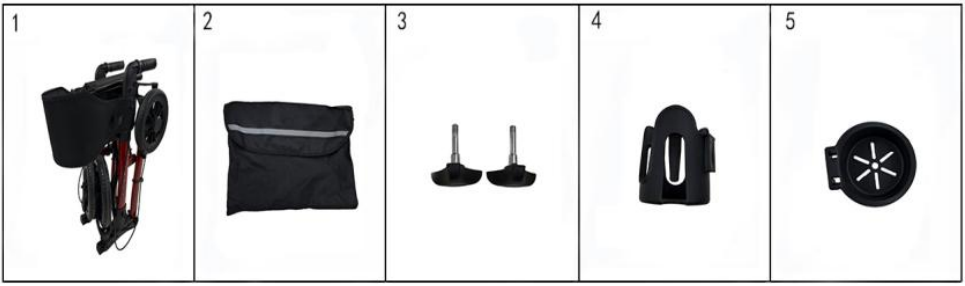
Figure 1

1. Main frame set
2. Hanging bag
3. 2 triangular knob screws
4. Cup holder
5. Cane box

Assembly Instructions: (Please refer to the part numbers and names in Figure 1 for assembly):

Step 1: 1. (As shown in Figure 2) Take out the 1. Main frame set, rotate the two front wheels to the bottom of the frame. Then take the 3. 2 triangular knob screws, align them with the holes on the front wheel assembly, insert the triangular knob screws and tighten them. Next, lift the pull pin on the rear wheel tube of the 1. Main frame set, and rotate the rear wheel tube until the pull pin locks into another limit hole.

Figure 2

Step 2: (As shown in Figure 3) Press down the seat cushion with both hands until a "click" sound is heard—this confirms the seat cushion is fully opened. Pull up the pull strap of the seat cushion to unlock and fold it.

Figure 3

Step 3: (As shown in Figure 4) Snap 4. Cup holder and 5. Cane box into the cup holder mounting base and cane box base on the frame respectively. Then take 2. Hanging bag and hang it on the plastic hook of the cross tube.

Figure 4

Step 4: (As shown in Figure 5) Press the spring-loaded ball on the side of the 1. Main frame set, then lift or press the armrest up or down until the spring-loaded ball locks into the next hole to adjust the armrest height.

Figure 5

Step 5: Pull up the handbrake to slow down or stop (as shown in Figure 6). Push the handbrake down until a "click" sound is heard, which locks the four-wheeled vehicle in place (as shown in Figure 7). To unlock the vehicle's self-locking function, simply pull the handbrake up.

Figure 6

Figure 7

Post-Installation Inspection:

1. Verify that the four-wheeled vehicle is assembled correctly by comparing it with the diagrams in the manual.
2. Confirm that all components of the four-wheeled vehicle are fully assembled, and check for any leftover parts.
3. Check if any part of the four-wheeled vehicle is loose; if so, tighten the screws and nuts.
4. Check if all external components of the four-wheeled vehicle are misaligned; make appropriate adjustments if needed.
5. Check if the brake functions properly, the folding operation works smoothly, and the steering moves freely.

SAFETY PRECAUTIONS

1. Do not let children play with the packaging unsupervised, as it poses a suffocation risk to children.
2. Keep this product away from children. It is not a toy—children may tip it over and injure themselves.

3. Position the four-wheeled vehicle as close to your body as possible, and keep your body as straight and forward-facing as you can.
4. The four-wheeled vehicle may only be used for its intended functions.
5. The maximum weight capacity of the four-wheeled vehicle is 158kg.
6. The four-wheeled vehicle may only be used when it is in good condition.
7. Before use, ensure the four-wheeled vehicle functions properly and all screws are tightened.
8. The seat of the four-wheeled vehicle is only for short-term rest and cannot be used as a replacement for a chair.
9. Do not stand on the four-wheeled vehicle, and do not use it to carry other people.
10. Exercise extreme caution when using it on road shoulders, near obstacles, on sloped, uneven, or slippery surfaces, as this may lead to falls or injuries.
11. The four-wheeled vehicle may only be used after a thorough inspection, and when it is fully deployed and locked in place.
12.Be careful not to get your hands caught when folding or deploying the vehicle.

PRODUCT MAINTENANCE AND CARE METHODS

1. Do not scratch or scuff the product surface with hard objects to avoid damaging the material.
2. For daily cleaning and maintenance, use a neutral synthetic detergent: wipe with a damp cloth first, then dry with a dry cloth.
3. Do not let the product come into contact with strong acids or alkalis to prevent corrosion.

STORAGE CONDITIONS AND METHODS

After packaging, the four-wheeled vehicle should be stored indoors in a clean, well-ventilated area with no corrosive gases, where the temperature stays within -10°C to 40°C and the relative humidity does not exceed 80%. During transportation, protect it from moisture, avoid exposure to light, and keep it away from heat sources.

CORRECT DISPOSAL

This product is subject to the provision of european Directive 2012/19/EU. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to acollection point for recycling electrical and electronic devices.

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