VEVOR

Upgrade · The Home Creator Way

# Tub Transfer Benches

## Model: FBL650329

VEVOR
Upgrade·The Home Creator Way

# Tub Transfer Benches

**MODEL:FBL650329**

<Picture Only For Reference >

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

# IMPORTANT SAFEGUARDS

**WARNING:** Read and understand this entire manual before operating or servicing this product. Failure to follow these warnings and instructions can cause personal injury or damage to valuable property.

# ASSEMBLY PRECAUTIONS

1. Assemble only according to these instructions. Improper assembly can create hazards. The manufacturer declines all responsibility in case of damage due to a bad installation.
2. Caution with the risks of pinching or cuts during assembly.
3. Avoid children using the product. And this product is not a toy. Do not allow children to play.
4. Wear ANSI-approved safety goggles and heavy-duty work gloves during assembly and use.
5. Keep the assembly area clean and well-lit.
6. Keep bystanders out of the area during assembly.
7. Do not assemble when tired or when under the influence of alcohol, drugs, or medication.
8. Product capabilities apply to a properly and completely assembled product only.
9. The product complies with relevant standards and regulations.Installation of this device must comply with all applicable regulations,and it should only be operated in well-ventilated areas.
10. This product is not intended for use by persons (including children) with reduced physical,sensory or mental capabilities or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
11. Please properly keep product. Check whether the product is abnormal before use. If it is abnormal, please stop using it.
12. Please use the product according to the intended use, and warn of potential injury from misuse.
13. For additional information regarding the parts listed in the following pages, please refer to the Assembly Diagram of this manual. Unwrap and separate all parts in a clean work area. Please keep small spare parts out of children's reach.

# SAVE THESE INSTRUCTIONS

## WARNING:



- Read the manual completely to avoid harm, product damage, or warranty avoidance. Seek professional help if you encounter discomfort, pain or worsening symptoms. If necessary, use supervision in order to assemble and use this product correctly and safely. Do not exceed the weight limit.
- Regularly check for secure assembly and contact support if in doubt. Ensure the product's legs are even and locked before use, and discontinue usage if there are signs of wear, damage, or missing parts. The product is to be used on flat, stable surfaces. Do not disassemble if not necessary, as it might compromise the longevity of the product. Do not modify the product unless specified in the manual.
- Exercise caution during operation and assembly. Avoid extreme temperatures or open flames. If an allergic reaction occurs, cease use immediately and consult a healthcare professional.
- Clean regularly using a mild detergent and warm water, rinse off any soap residue. After cleaning, dry thoroughly
- Not suitable for environments with high humidity, explosive atmospheres, salt spray, acidic or alkaline conditions.

## IMPORTANT


1) Don't exceed the maximum weight capacity.
2) Only for use on bathtubs that are stable and have flat surfaces on both sides.
3) Regularly inspect for signs of wear, damage, or loose parts.
4) Before use, confirm with your healthcare provider that the product meets your specific needs. The following guidelines are not a substitute for medical advice.

## PART LIST

When unpacking this product, check to make sure the parts listed below are included and carefully inspect for any damage that may have occurred during transit. If any part is missing or damaged, do not attempt to assemble or use the product.
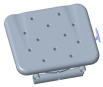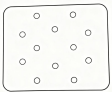
| NO | PART NUMBER | PICTURE | QTY |
|---|---|---|---|
| **1** | Backrest (with safety belt) |  | 1 |
| **2** | Seat board set |  | 1 |
| **3** | Seat frame tube |  | 2 |
| **4** | Trestle A |  | 1 |
| **5** | Trestle B |  | 1 |
| **6** | Chair foot |  | 4 |
| **7** | Armrest |  | 1 |
| **8** | Shower buckle |  | 1 |
| **9** | Seat pad |  | 1 |
| **10** | Backrest pad |  | 1 |
| **11** | Manual |  | 1 |

# DIMENSIONAL PARAMETERS
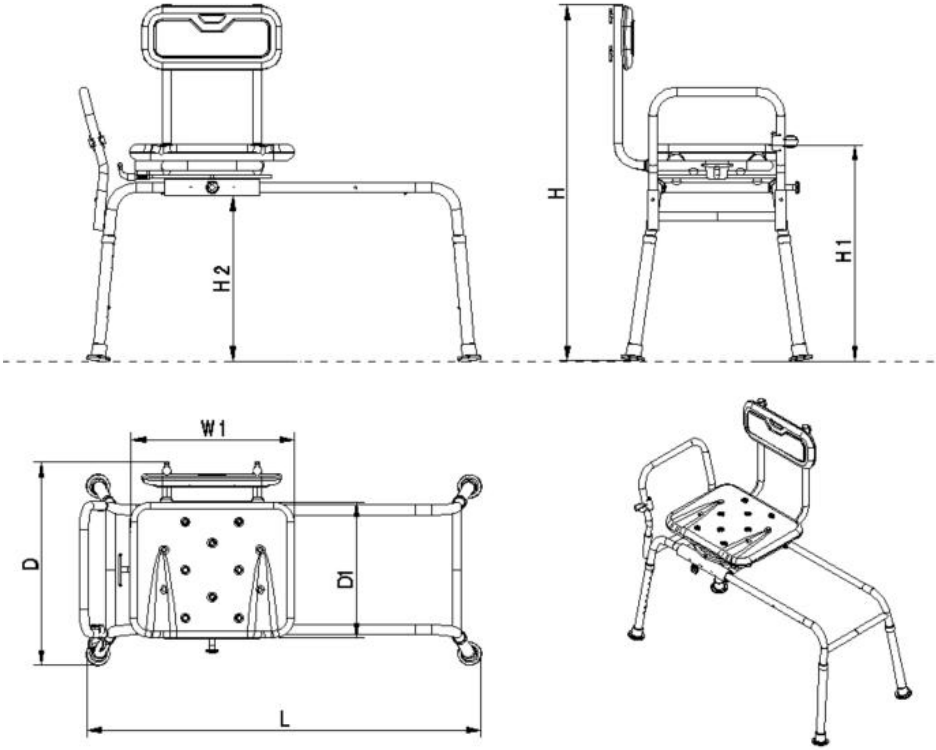
| Unit:(mm/in) | | | |
|---|---|---|---|
| L | D | H | D1 |
| 1000-1035<br>(39.4-40.7) | 520-540<br>(20.5-21.3) | 845-1020<br>(33.3-40.2) | 340(13.4) |
| W1 | H1 | H2 | |
| 410<br>(16.1) | 525-700<br>(20.7-27.6) | 380-555<br>(15-21.9) | |


**※There may be a ±20mm error in the measurement of some dimensions, which is normal and not a product quality issue.**

## PAREMETERS


| Model | FBL650329 |
|---|---|
| Main Materials | Aluminum |
| Rotation Angle | 360° |
| Seat Rated Load | 180kg/400lbs |


# INSTALL INSTRUCTIONS

**Step1**.Thread the two Seat frame tubes through the two tubes on the front and rear sides of the Seat board set in sequence
(※**Press the button while assembling the tube**).
**Push the tube until you hear a 'click' sound, indicating that it is securely locked.**

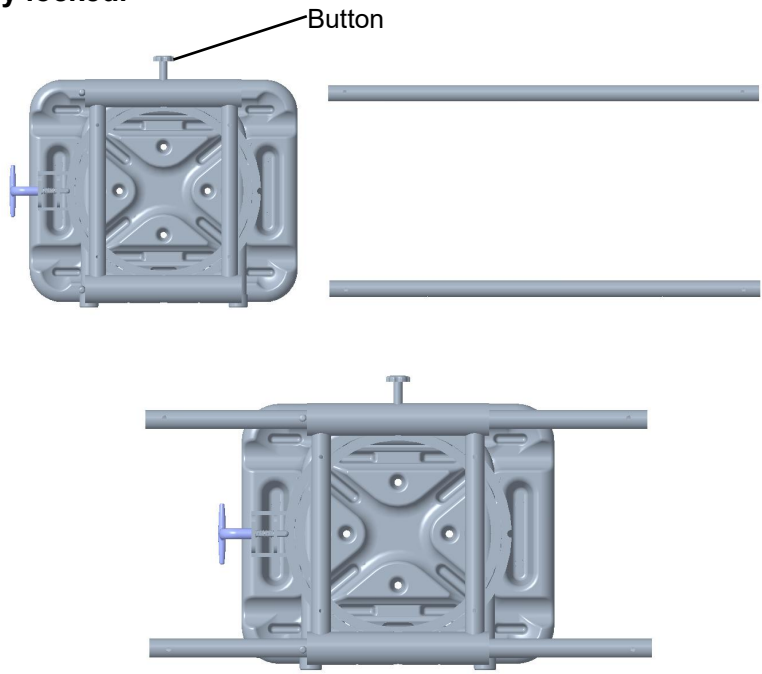
**Step2**.Press the ball pin on the leg tube and insert it into the hole on the side of the trestle(A&B),Repeat on both sides.

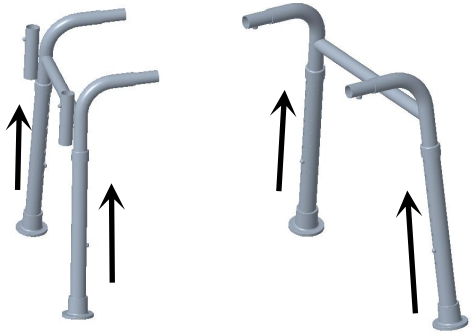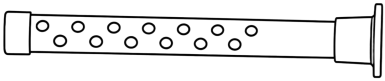The chair legs are adjustable in 15 gears, and users can adjust the position of the chair leg ball hole to meet the height requirements of the seat.

**Step3**.Align the legs(part of step2) insert them into the Seat frame tube (ensure that the holes of the bracket pipe openings are aligned with the holes of the bracket pipe openings).

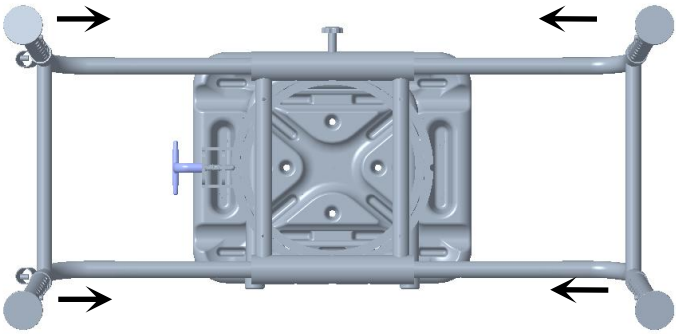**Flip it over and make sure the four foot pads are standing on the same level.**


**Step4**.Align the backrest with the pipe opening on the back of the seat frame , and then insert it. (ensure that the ball pin of each backrest pipe pops out normally)


**Step5**.Connect the armrest to the bracket on the side.

**Step6**.Fix the shower buckle to the armrest.

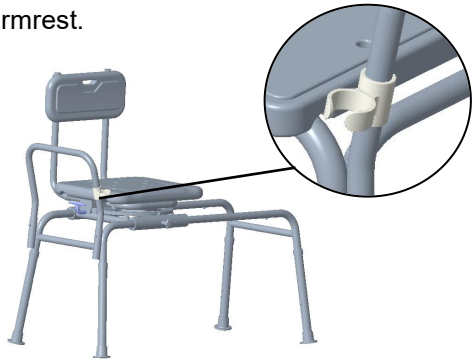To adjust the height, press the ball pin on each leg and slide the leg up or down to your desired height. Ensure all legs are set to the same height and lock securely into place.

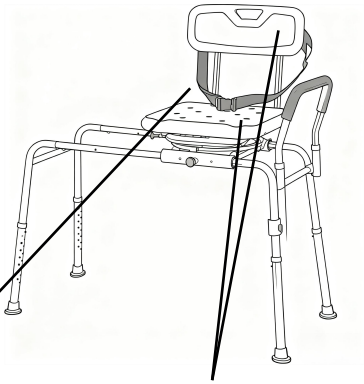
The seat belts was fixed each of the two backrest to secure the user's body.

Attach the included backrest pad and seat pad for a more comfortable for use.

## USE INSTRUCTIONS


## 1.Leveling

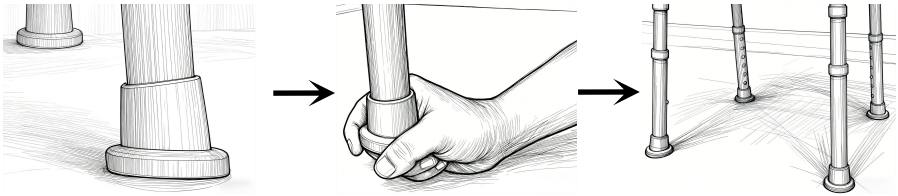If the foot mat is not placed flat on the ground, rotate it until it is fully in contact with the ground.


## 2.Swivel the Seat.

After the product assembly is completed, a pull rod can be seen below the left side of the seat plate. When pulling the pull rod, the direction of the seat plate can be adjusted by rotation. After releasing the pull rod, the seat plate will rotate in one direction and fix its direction by itself.
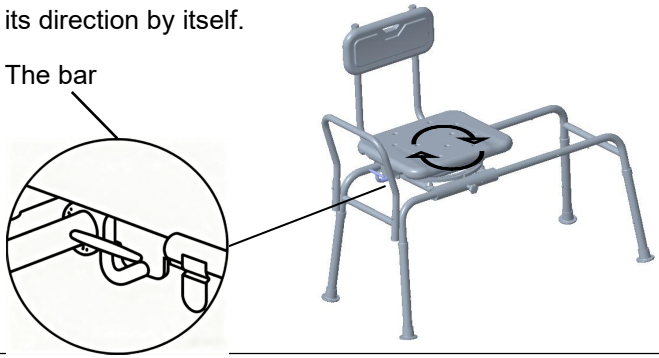
## 3.Slide the seat.

Pull out the button under the seat board, and the seat can slide left and right. When the seat is aligned with the hole on the main frame tube, you will hear a "click" sound. At this point, the seat board is firmly locked and can be used normally.

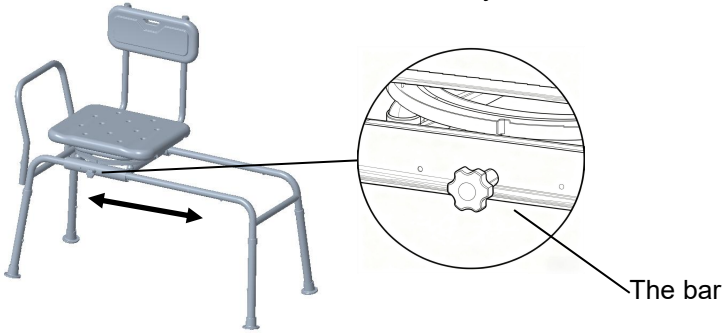**USING WARNING:**

(1) Put the shower bench stably on a flat floor.

(2) Before use, please make sure that the screw has been tightened and the ball pin of the legs has been ejected from the slot.

(3) Before use, adjust the shower bench to an appropriate height so that when the two legs are supported in the bathtub, the seat is smooth.

(4) Check the stability of the shower bench before using it every time. Stop using it if the legs are shaking or the shower bench is unstable.

(5) Avoid the overload.

(6) The shower bench should not be used to push the user while seated.

(7) Shower bench must be used ONLY for taking a shower inside bathroom, not for any other purposes.

(8) Adjust the shower bench legs to a suitable height so that the user's feet from toes to heels comfortably touch the floor when sitting on it.

(9) Stay away from open flames.

(10) Periodically check whether screws have been tightened.

(11) Ensure that rubber tips at the bottom of the shower bench’s legs are not worn out. If not, replace them immediately.

(12) Please read carefully before use and do not tear off the label.

# MAINTENANCE

- Wash with clean or mild soapy water and air dry after use.
- When not in use, please store in a cool and dry place
- Check the shower bench regularly.

# STORAGE CONDITIONS

* The fabric should be stored completely dry, in a clean and dry place to avoid the risks of mold.

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